

TECHNICAL SPECIFICATION*

Throughput	: 200 photometric tests per hour and 400 tests per hour with ISE* (Optional ISE with Na ⁺ , K ⁺ , Cl ⁻ , Li ⁺)
System Type	: Discrete, Open, Automated, Random Access, Patient Prioritized Clinical Chemistry Analyzer.
Onboard Parameters	: upto 50 + 4 ISE*
Sample Type	: Serum, Plasma, Urine, CSF, HbA1c using Whole Blood*
Programmable Parameters	: Unlimited
Analytical Methods	: 1-Point, 2-Point, Rate-A, Rate-B, Direct Potentiometry
Calibration	: Linear, Non-Linear, Multipoint
Photometer	: Static Photometer
Absorbance Range	: 0-3.0 Absorbance
Light Source	: Halogen Lamp 12 V/ 20 W
Optics	: 8 Filters (340-700nm) 340, 405, 505, 546, 578, 600, 660, 700nm.
Detector	: 8 Silicon photo diodes
Sample Unit	: 39 positions for samples/ blank/ calibrator / control & STAT sample.
Reagent Unit	: 50 cooled reagent positions
Quality Control	: QC Plot and data with QC rules
Reaction Tray	: 45 Hard glass Cuvettes
Reaction Liquid Mixing	: Stirrer with variable speed Mixer
Reading Volume	: 180 µl.
On board laundry	: 6 stage cleaning, 2 stage Drying with Cuvette validation step
Sample Pipetting	: 2-70 µl (Adjustable in 0.1 µl)
Reagent Pipetting	: 10-300 µl (adjustable in 1 µl)
PC Configuration	: OS-Windows 7 embedded or higher, CPU -P4 or higher, RAM 2GB, HDD - 80GB, USB connectivity
Power supply	: AC 110 V.+/- 10% 60±1 Hz or AC 220 V ± 10%,50± Hz.(Factory Set)/600VA
Water Consumption	: Upto 7.5 ltrs per hour
Hibernate	: Enhances Lamp life and Pump life
Dimensions	: 810 mm(W) x 800 mm(D) x 600 mm(H)
Weight	: 130 Kg

* Using onboard lyse features supported by XL 200

* Specifications subject to change without prior notice

XL200

Automatic Clinical Chemistry Analyzer



Optimal solution
for small and medium
laboratories



PRM Medical (Pty) Ltd
76 Skilpad Avenue, Monument Park, Pretoria, 0181
087 265 7919
info@prmmmedical.co.za

PRM
Medical
Clinical Diagnostic Solutions Advisory and Supplier

XL200

Available automation of clinical chemistry analysis

DISPENSING OF SAMPLES AND REAGENTS

- Sample volume: 2-70 µl (in 0,1 µl step)
- Reagent volume: R1 50-300 µl (in 1 µl step),
R2 10-200 µl (in 1 µl step)
- Dispensing probe equipped with liquid -level sensor and crash detector
- Auto-dilution of samples and calibrators

ECONOMY



- Minimum reaction volume: 180 µl
- Reusable reaction cuvettes

MIXING SYSTEM

- Independent stirrer
- 3 user selectable mixing speeds

QUALITY CONTROL

- 4 levels of control material can be used
- Levey-Jennings graphs
- Twin Plot diagrams for monitoring of systematic and random error

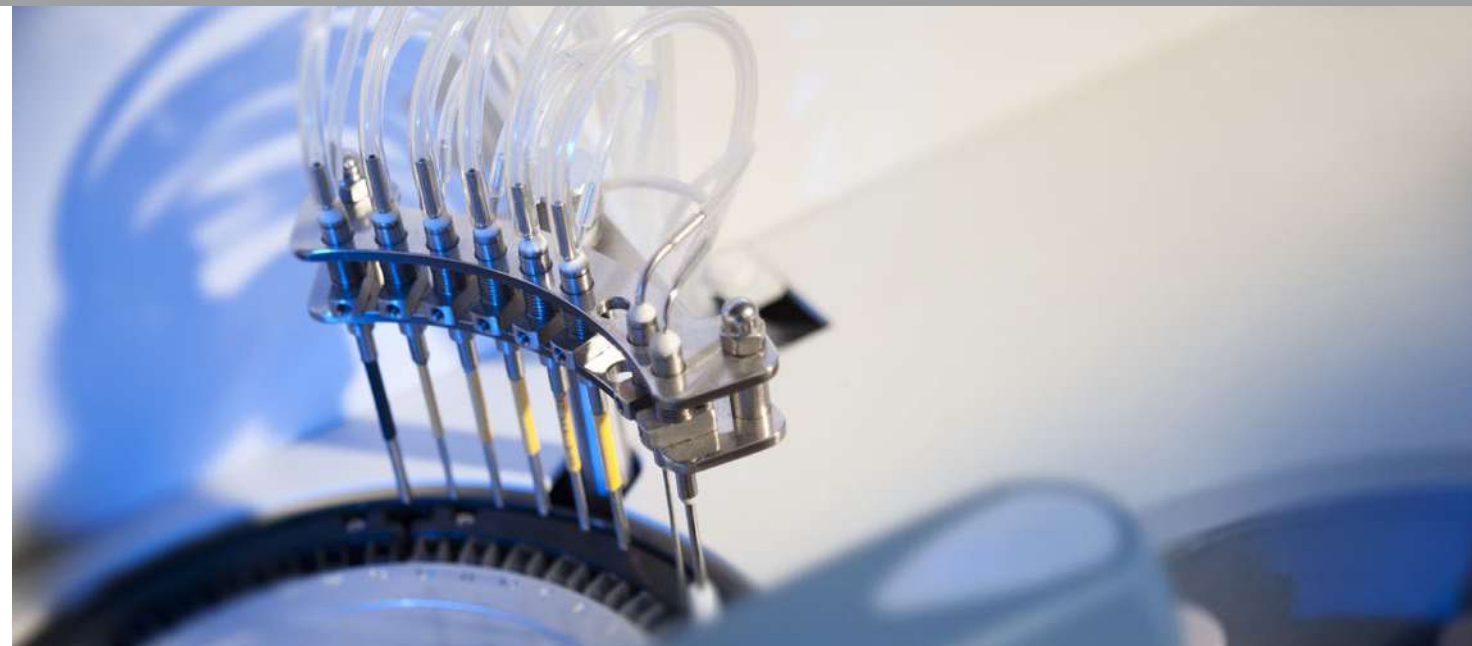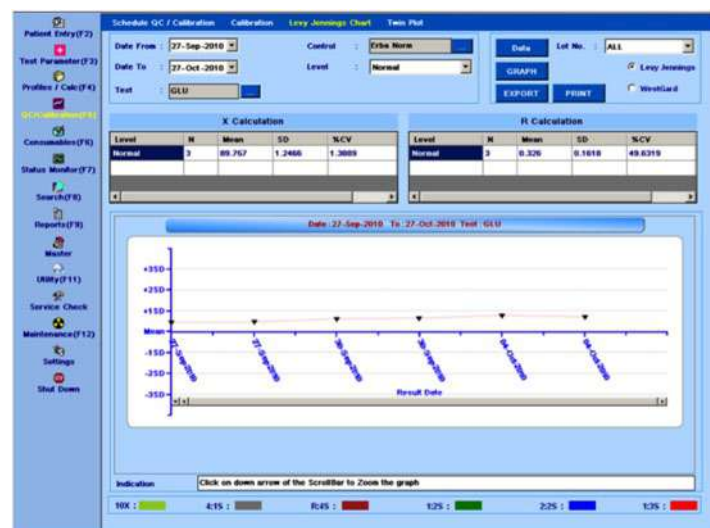

REACTION UNIT

- 45 reusable hard glass cuvettes
- Possibility of replacement of individual cuvette
- Wash station cuvette rinsing and drying in 8 steps
- Automatic cuvette blank measurement before analysis
- Reaction temperature 37°C ± 0,2°C

SAMPLE TRAY

- 39 positions for samples, blanks, standards, calibrators, controls and ISE solutions
- Primary tubes 5, 7 and 10 ml and cups
- STAT sample with priority in any position
- Possibility of programming up to 99 virtual trays



REAGENT TRAY

- 50 positions, 20 ml, 50 ml reagent containers, 5 ml tube (with adaptor)
- Reagent compartment with Peltier/air cooler (8-12°C)
- Option to use one reagent for several tests simultaneously

SOFTWARE

- Convenient user interface
- Connection to LIS
- Statistical methods of processing results
- Data export in selected format



MEASUREMENT MONITORING

- Color indication of sample analysis
- Possibility of monitoring the reaction in real time
- Reagent volume monitoring
- Informative reports on ongoing analyzer status

